

LED Phototherapy Unit

JW-PU1000

CE

2460

Presafe

ALONE & REMOVED

COMPONENT



LED Phototherapy Unit


JW-PU1000

Functionality

- 3 steps brightness function
- Timer function
- Operating time and total usage time display


Mobility

- Portable and compatible with JW infant care product series

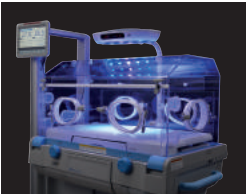


Effectiveness

- Consisting 2 white LEDs and 8 blue LEDs
- Blue LEDs -->Effective range of light treatment (400~550nm)
- White LEDs -->Examination purpose

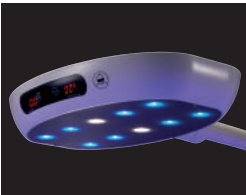

Adjustability

- Wide exposure range by tilting and height adjustment

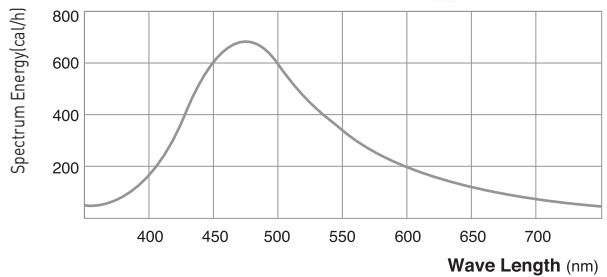


Blue LEDs emit light in the range of 400 – 550 nm is considered to be the most effective for the degradation of bilirubin

Specification

CONTENTS	DESCRIPTION
Dimension(mm)	Housing: 377 x 274 Base: 403 x 428
Height Adjustment(mm)	1200~ 1700
Lamp	2 White LED, 8 Blue LED
Power Source	AC 100~240V, 50/60Hz
Power Consumption	45VA
Set Time	30~120 minutes (interval of 30 minutes)
Spectral Irradiance	60uW/cm²*nm
Caster	65mm rotary wheel x 4EA
Housing Tilt	±35°
Weight	24kg
Lamp Life	20,000 hours

Standard Accessory

CONTENTS	EA
Eye Mask	Large 1, Small 1
Curtain	1
Manual	1
Power Cord	1

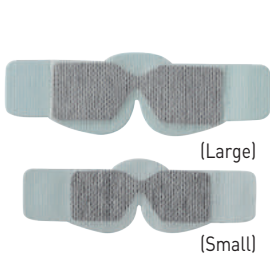

Eye Mask



Curtain